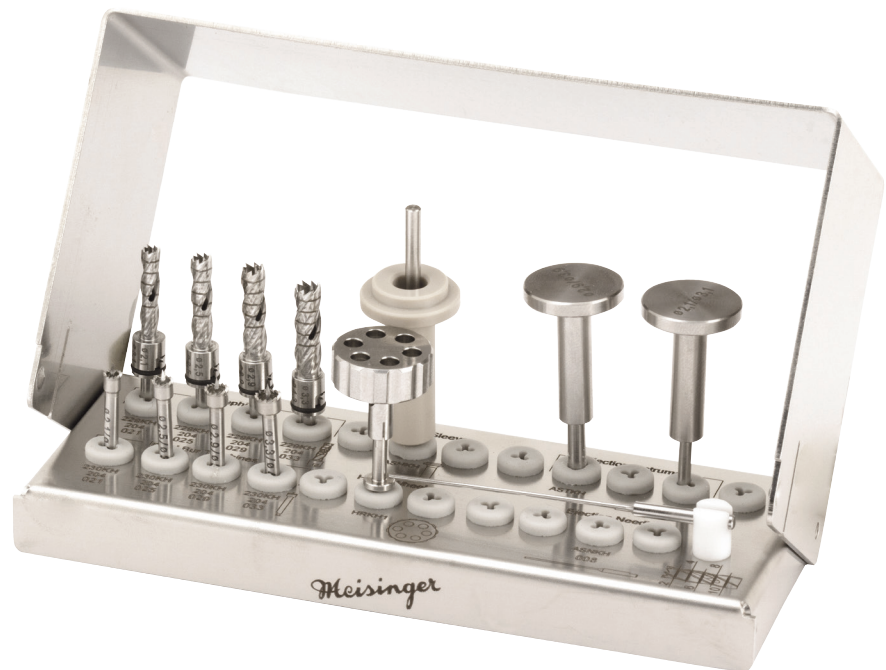


Trephine Ejection Kit

By Prof. Dr. Fouad Khoury

The Trephine Ejection Kit combines several unique ideas, which makes this set something very special. Internally cooled, the trephines provide a safe extraction of bone cylinders. The internal cooling protects the bone tissue from damage due to overheating. The removable trephine working parts allow for a particularly easy extraction of the cut

out bone cylinders. With the aid of special ejection instruments, which fit to the trephine diameters, the bone can be pushed out easily. Further, the kit contains initial trephines which provide the safe application of the trephines. The Trephine Ejection Kit allows for fast and safe working and provides high flexibility due to various instrument sizes.



BTE00

At a glance



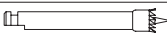
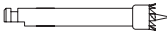
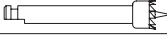
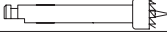
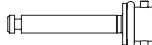




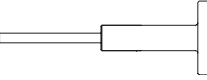
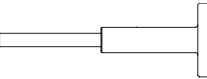
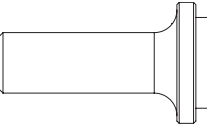
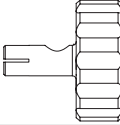

- Two-part trephines for a safe and minimal invasive extraction of bone cylinders
- Protection against overheating due to internally cooled trephines
- Innovative application
- Standardized diameters and precise laser markings of the trephines allow for a predictable extraction of defined bone cylinders
- Perfectly matched instruments care for highest precision and flexibility

Our sales force would be pleased to answer your questions: +49 2131 2012-303

Meisinger since 1888

Trephine Ejection Kit

Fax for ordering +49 2131 2012-222

| | Figure | Shank | Size | Internal / External diameter | recomm. / max. revolutions | Pieces / Kit | Amount  |
|---|--------|-----------|------|------------------------------|---|--------------|--|
|  | BTE00 | - | - | - | - | - | |
|  | 230KH | RA (204) | 021 | 2.1 mm / 3.1 mm | 800 min ⁻¹ / 800 min ⁻¹ | 1 | |
|  | 230KH | RA (204) | 025 | 2.5 mm / 3.5 mm | 800 min ⁻¹ / 800 min ⁻¹ | 1 | |
|  | 230KH | RA (204) | 029 | 2.9 mm / 3.9 mm | 800 min ⁻¹ / 800 min ⁻¹ | 1 | |
|  | 230KH | RA (204) | 033 | 3.3 mm / 4.3 mm | 800 min ⁻¹ / 800 min ⁻¹ | 1 | |
|  | 229KH | RA (204) | - | - | 400 min ⁻¹ / 4.000 min ⁻¹ | 4 | |
|  | 229KH | unmounted | 021 | 2.1 mm / 3.1 mm | 400 min ⁻¹ / 4.000 min ⁻¹ | 1 | |
|  | 229KH | unmounted | 025 | 2.5 mm / 3.5 mm | 400 min ⁻¹ / 4.000 min ⁻¹ | 1 | |
|  | 229KH | unmounted | 029 | 2.9 mm / 3.9 mm | 400 min ⁻¹ / 4.000 min ⁻¹ | 1 | |
|  | 229KH | unmounted | 033 | 3.3 mm / 4.3 mm | 400 min ⁻¹ / 4.000 min ⁻¹ | 1 | |
|  | ASTKH | - | 021 | - | - | 1 | |
|  | ASTKH | - | 029 | - | - | 1 | |
|  | ASHKH | - | - | - | - | 1 | |
|  | HRKH1 | - | - | - | - | 1 | |
|  | ASNKH | - | 008 | - | - | 1 | |

CE 0044

Available from the dental trade.

Your trade depot:

Your delivery address:

Hager & Meisinger GmbH

Hansemannstr. 10 • 41468 Neuss • Germany
 Phone: +49 2131 2012-0 • Fax: +49 2131 2012-222
 info@meisinger.de • www.meisinger.de

Meisinger USA, L.L.C.

10200 E. Easter Avenue • Centennial, Colorado 80112 • USA
 Phone: +1 (303) 268-5400 • Fax: +1 (303) 268-5407
 Toll free: +1 (866) 634-7464
 info@meisingerusa.com • www.meisingerusa.com

85FL119EN - 0120