

Occlusal Reduction Kit & Tapered Occlusal Reduction Kit

The Occlusal Reduction Kit allows precise occlusal reductions to be carried out. Thanks to the unique, angled depth stops a precisely defined removal is possible. The vertical orientation provides for good access to fissures. A crosshatched pattern is formed on the occlusal surface that, when connected, yields a perfect inclined plane reduction. The diamond burs are available in 1.5 mm, 1.8 mm, 2.0 mm, 2.2 mm, and 2.4 mm length corresponding to proper occlusal reduction for various restorative materials. A 1.0 mm bur is also available for labial reduction. The Occlusal Reduction Kit contains cylindrical diamonds while the instruments of the Tapered Occlusal Reduction Kit are conical.

Made in Germany



2565



2565T

At a glance










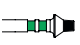

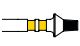
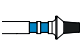
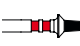
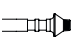
- Precise occlusal reduction
- High defined removal
- Easy handling
- Depth variations for adaptation to various restoration materials

Our sales force would be pleased to answer your questions: +49 2131 2012-303

Meisinger since 1888

Occlusal Reduction Kit & Tapered Occlusal Reduction Kit

Fax for ordering +49 2131 2012-222

	Figure	Shank	Size	recomm. / max. revolutions	Packaged à	Amount 
	2565	-	-	-	-	
	828G	FG (314)	017	160.000-250.000 min⁻¹ / 250.000 min⁻¹	1	
	828O	FG (314)	017	160.000-250.000 min⁻¹ / 250.000 min⁻¹	1	
	828Y	FG (314)	017	160.000-250.000 min⁻¹ / 250.000 min⁻¹	1	
	828B	FG (314)	017	160.000-250.000 min⁻¹ / 250.000 min⁻¹	1	
	828R	FG (314)	017	160.000-250.000 min⁻¹ / 250.000 min⁻¹	1	
	828W	FG L (315)	017	160.000-250.000 min⁻¹ / 250.000 min⁻¹	1	
	2565T	-	-	-	-	
	828TG	FG (314)	017	160.000-250.000 min⁻¹ / 250.000 min⁻¹	1	
	828TO	FG (314)	017	160.000-250.000 min⁻¹ / 250.000 min⁻¹	1	
	828TY	FG (314)	017	160.000-250.000 min⁻¹ / 250.000 min⁻¹	1	
	828TB	FG (314)	017	160.000-250.000 min⁻¹ / 250.000 min⁻¹	1	
	828TR	FG (314)	017	160.000-250.000 min⁻¹ / 250.000 min⁻¹	1	
	828TW	FG (314)	017	160.000-250.000 min⁻¹ / 250.000 min⁻¹	1	

CE 0044

Available from the dental trade.

Your trade depot: <hr/> <hr/> <hr/>	Your delivery address: <hr/> <hr/> <hr/>
--------------------------------------------------	-------------------------------------------------------

Hager & Meisinger GmbH
 Hansemannstr. 10 • 41468 Neuss • Germany
 Phone: +49 2131 2012-0 • Fax: +49 2131 2012-222
 info@meisinger.de • www.meisinger.de

Meisinger USA. L.L.C.
 10200 E. Easter Avenue • Centennial, Colorado 80112 • USA
 Phone: +1 (303) 268-5400 • Fax: +1 (303) 268-5407
 Toll free: +1 (866) 634-7464
 info@meisingerusa.com • www.meisingerusa.com

Please observe the general instructions for the use and safety of MEISINGER products in the medical field as well as the notes on reprocessing (cleaning, disinfection and sterilization) Hager & Meisinger GmbH medical devices.

85FL115EN - 0317